

Treating hands like flowers



External Fixator
with dorsal approach



EXTERNAL
fixation



monolateral
fixators



hand

Indications

Indicated for gradual distraction of inter-phalangeal and metacarpo-phalangeal joints in:

- post-traumatic joint contractures
- high degree Dupuytren contractures






MiniFlo R25



MiniFlo R30



Offset Extension

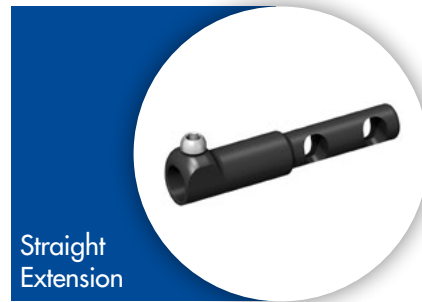

Straight Extension

Features	Benefits
Arc shaped fixator body	Easy dorsal application on the hand
Device positioning over the rotation axis of the joint	Gradual symmetrical distraction without stressing joint structures
Dual function Targeting device	Convenient device for correct fixation positioning intraoperatively and as a distraction wrench postoperatively
Easy distraction mechanism	Patient-friendly: the patient can easily distract on the basis of tissue tension (as tolerated)
Sterile Bone screws available in three sizes	Versatility of use for different bone dimensions
Fixator available in two sizes	Application on inter-phalangeal and metacarpo-phalangeal joints
Supplementary Telescopic Extensions	The extension prevent asymmetric or over-distraction when two joints are involved, or in case of over distraction

CODE	DESCRIPTION	STERILE	
F4R-1025	MiniFlo Fixator R25 (supplied with MiniFlo Key, 1.5 mm K-wire and 3 mm Allen wrench)		
F4R-1030	MiniFlo Fixator R30 (supplied with MiniFlo Key, 1.5 mm K-wire and 3 mm Allen wrench)		
F4R-2021	MiniFlo Offset Extension (only for F4R-1030)		
F4R-2022	MiniFlo Straight Extension (only for F4R-1030)		

Bone Screw, total length 100 mm, thread length 13 mm		STERILE	
SF32160	Stem \varnothing 2.5 mm, thread \varnothing 1.6 mm		
SF32200	Stem \varnothing 2.5 mm, thread \varnothing 2.0 mm		
SF32250	Stem \varnothing 2.5 mm, thread \varnothing 2.5 mm		